seca incscale 748

## Continuous operation

When the mains adapter unit is connected, the scale will automatically operate continuously for 10 minutes.

## Tare function

It is possible to tare objects such as blankets and covers before weighing:

- Place objects on the scale.
- Press Tare button (8) – "0.000" will appear and the tare indicator (9) above the button will light up.
- Actual weighing can now take place. The tare function can be used as often as you like up to maximum load. The maximum weight indicated will be reduced by the objects already tared.
- If an object is removed from the scale during the tare procedure, then "-. - -" will appear. Press the tare button (8) and "0.000" will reappear.
- When there is absolutely no load on the scale, the tare function stops automatically. "0.000" will be indicated and the tare indicator light goes out.

## What do I do if...

... "- - - -" appears in the indicator?

- Press tare button (8).
- Restart scale (10).

... "5E OF" appears in the indicator?

- The scale is overloaded.

... "EEEE" appears in the indicator?

- The automatic monitoring system has discovered a malfunction somewhere

in the electronics.

... "bAtt" appears in the indicator?

- The battery needs to be recharged or replaced (11) (Recharging time approximately 24 hours).
- Operate the scale with the mains adapter unit.

## Function security check of unloaded scales

Measurement function security of the unloaded scale can be checked by pressing the Check button (7). The indicator will show the word "tEt". After this, the luminous displays and all the segments of the weight indicator "0000" will flash on and off several times. The final function security figure to be indicated must lie within the range 8995 to 9005. The word "End" indicates that the test is completed.

## Technical Data

Weighing range:	0–9.995 kg
Tare range:	–9.995 kg
Fine scale:	5 g
Maximum load permitted (baby + cot bed):	15 kg
Manual tare function	
Function check via test values	
Automatic readjustment to zero	
Stable indications due to electronic damping	
Automatic reset to zero at start after the baby has been lifted	

Power supply:

Battery operation for about 60 weighings  
9 V NiCd square battery

Tr 7/8 9 V 100 mAh

Mains operation with simultaneous battery recharging Adaptor

220 V = 50 Hz

11 V = 70 mA

68-32-10-252

seca – Order No.:

Dimensions:

Height 210 mm

Width 850 mm

Depth 350 mm

Weight: approx. 5 kg

Calibration: medical calibration class III

## What else you need to know...

- Use the customary disinfectants and cleaning agents to clean the scale.

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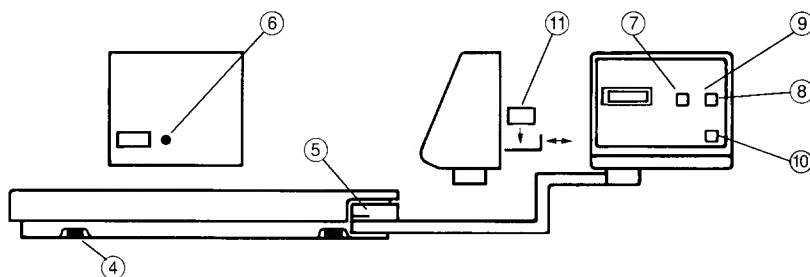
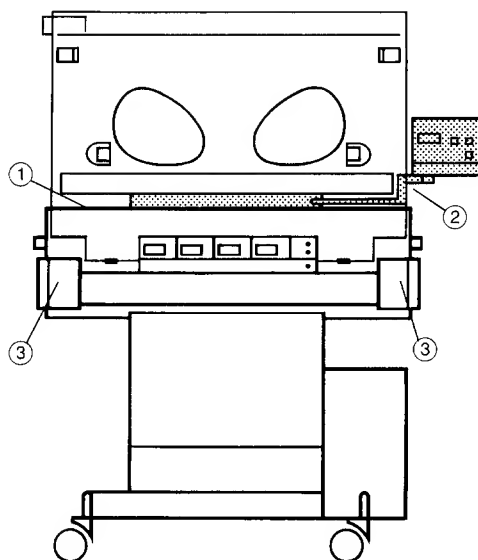
## 1 year

warranty is given after date of purchase for all faults due to material or construction (except for batteries and electrodes). Faults under this warranty are repaired without charge. Any further claims will be excluded. Transport both ways will be to the charge of the customer if the equipments are not at the place of the purchaser. Warranty claims will be invalidated if the equipments have been tampered with by unauthorised personnel.

Customers abroad are asked to contact the seller in their country for warranty claims.

**Important:** Warranty can be claimed only together with the purchase voucher.

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